


# P312A [DUALER]

Convertible FWD/RWD

Powerbase Wheelchair

 **merits**  
**AUSTRALIA**

[www.meritsaustralia.com](http://www.meritsaustralia.com)

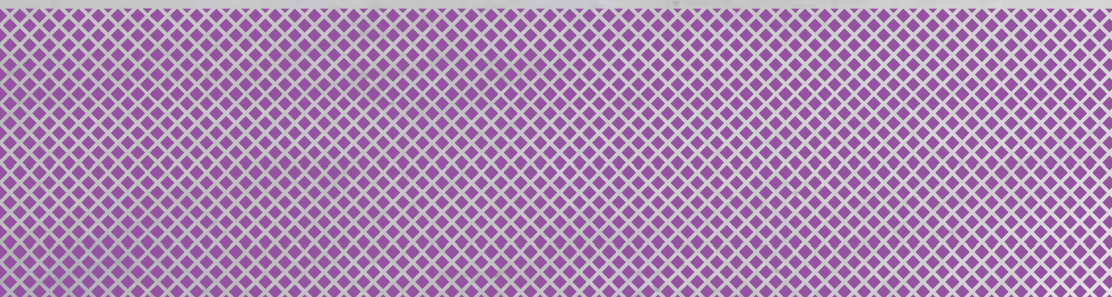



## Features

### Unique and compact

- 2-in-1 powerchair, easy to convert from RWD to FWD, ideal chair for indoor as well as outdoor use.
- Unique compact size for maximum manoeuvrability inside the home or in narrow spaces and walkways.
- Height & width adjustable flip-up armrests for improved comfort.

### Stable and flexible

- Change from RWD and move to FWD by simply rotating the chair 180°, to have the perfect drive option for every situation.
- Turning radius is only 49cm/ 19.6”.
- Low seat height of 18” provides excellent stability.
- Power seat elevator available.



Power seat elevator available.  
13cm/5in



Very easy to convert from rear to front wheel drive: Just release and rotate the chair 180°

## Dimension (cm/in)

Seat Width	Seat Depth	Seat Height from Ground	Seat Height from Deck	Backrest Height	Overall Length Footplate Down	Overall Width	Overall Height
46/18	46/18	56/22	46/18	51/20	101/40	61/24	104/41

## Technical Specs

Weight Capacity	Max.Speed	Range	Turn Radius	Ground Clearance	Motor	Controller	Battery	Charger
136kg/300lb	8kph/5mph	28km/18mi.	53cm/21"	5cm/2"	DC 24V, 180W(Rating)/ 630W(Max)	PG VR2 60A	12V/ 34Ah x 2 pcs	5A off-board
Gradient	Caster	Drive Wheel	Brake	Base (kg/lb)	Battery/ea (kg/lb)	Seat (kg/lb)	tll w/o Batt. (kg/lb)	tll w/Batt. (kg/lb)
12°	8"(200x50) PU tire	10"(260x85) foam filled tire	Intelligent, regenerative, electromagnetic brakes	44/97	12/26	18/40	62/137	86/189